


URETHRAL CATHETERISATION DEVICE (UCD®)



Demo Video

Scan QR Code

or visit <https://scvn.io/dV0C>

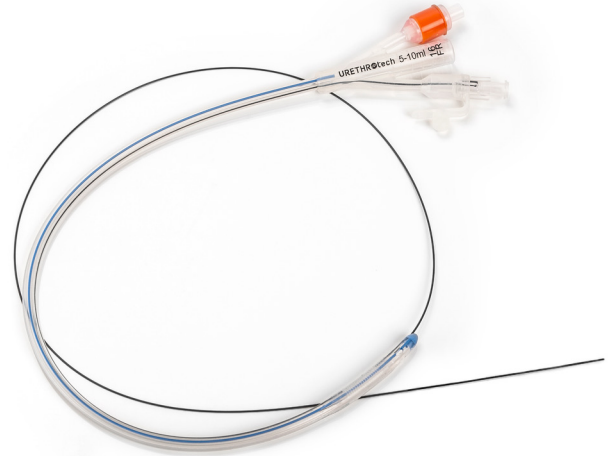

Device Description

The innovative Urethrotech UCD® (Urethral Catheterisation Device) is composed of a 3-way Silicone Foley Catheter with integrated introducer guidewire for safe bladder access to prevent Catheterization-associated Urethral Injury (CAUI). The integrated removable guidewire is made of Nitinol with hydrophilic coating to facilitate the atraumatic introduction of the catheter through the urethra into the bladder to negotiate any prostatic curve or tight external urethral sphincter. Also suitable for irrigation of the bladder after guidewire removal.

- > Completely 100% Silicone 3-way Foley catheter
- > Clear for easy visualization of urine flow or debris
- > Radiopaque blue tip and strip for easy visualization during scans
- > Two large opposing eyes for efficient urine drainage
- > Integrated atraumatic hydrophilic Nitinol guidewire
- > Guidewire exits at tip of catheter for ideal catheter guidance
- > Secure guidewire Luer-lock to flush the guidewire with sterile water
- > Attached plug ready to seal guidewire channel after guidewire removal
- > Symmetrical balloon shape for added reliability

Urethral Catheterisation Device (UCD®)

3-Way 100% Silicone Foley Catheter with integrated guidewire

PACKAGING 1box/each

Item No.	Size	Balloon	Effective Length	HCPCS/CPT	Pkg.
UCD-0035-090-0016	16 Fr	10 mL	90cm	A4346, 51703	1/bx



Specifications

HPIS CODE	730_10_40_20
BALLOON SIZE (CC)	10
CATHETER MATERIAL	Silicone, 100% Silicone
COATING	No
FRENCH SIZE (FR)	16
PRODUCT STYLE	3-way Foley Catheter, UCD
INTEGRATED GUIDEWIRE (CM)	90
UNSPSC	42142702

Unit of Measure	Conversion	Net Weight (g)	Volume (cubic m ³)	Shipping Dimensions L x W x H (cm)	GTIN-14
Box (BX)	1.0 Each	145	0.014742	52.5 x 10.8 x 2.6	15060639870008
Each (EA)	1.0 Each	145	0.014742	52.5 x 10.8 x 2.6	15060639870008

U.S. Sales Office:

Urethrotech Inc.
111 N Chestnut Street, Suite 104
Winston-Salem
North Carolina 27101-4054
Phone: 800-474-0180
Fax: 336-447-6085
orders@urethrotech.com

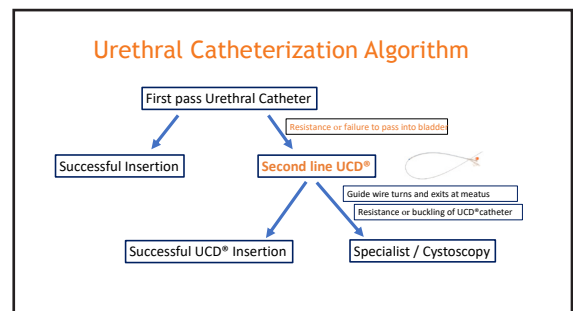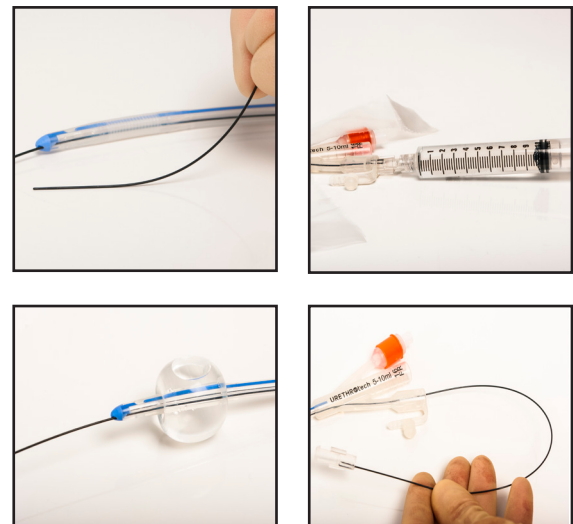

Urethrotech Ltd
24 Garth Road, Kingston Upon Thames, Surrey KT2 5NY, UK
www.urethrotech.com